


Single use and available sterile in boxes of 50 it can be used within the Operating Room at the point of the procedure to remove the full range of surgical blade shapes from either a No.3 or a No.4 fitment handle.

Protecting staff from sharps injuries during contaminated blade removal it is also transparent allowing for the count at the end of the procedure before disposing in an appropriate sharps container.

In stainless steel, with 5cm long cutting edge. Fits a No.3 fitment handle. Sterilised, packed individually, boxed in 20's.

Can be used as a component part of a "Procedure Kit". Small polystyrene handle can be supplied fitted with blade No's 10, 10A, 11, 15 or 15A. Polyethylene blade guard. Non-sterile and bulk-packed.

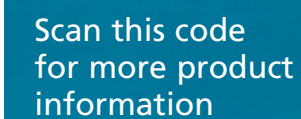

A choice of 3 'Sabre' sterile disposable scalpels – E/11, D/15 and B/23 – each supplied individually packed, boxed in 10's.

Available with the most popular stainless steel blade shapes, and are fitted with a polyethylene protective guard. The sterile scalpels are individually wrapped in a view pack, and boxed in units of 10. Also available in a bulk, non-sterile format. For use in custom procedure trays and kits.

Available with all standard stainless steel blade shapes, complete with CE marked polystyrene handles with a graduated measure (for indication only), individually sealed in a tray designed to assist with safe passing protocols, sterilized and boxed in units of 10.

Specially designed to assist in the reduction of potential sharps injuries this scalpel has a high visibility red handle combined with a retractable blade to support safe passing protocols. Colour coded sliders to indicate the blade shape, an audible click when engaging the usage, passing and locked positions and ribbed, slip resistant finger grips for a secure hold in a gloved hand. They are available in stainless steel blade shapes 10, 11P, 15, 20, 21, 22, 23, 24 and 25A in sterile packs of 25. Also available in non-sterile packs of 500. For use in custom procedure trays and kits.

The 'Fine' range has small, slim, precise blades in stainless steel for delicate operations. Six special handles, including the SF13 and SF23 with independently screwed top sections to facilitate blade-changing, take all blades, supplied in individual metal foil packets, available in dispenser units of 25. Handles are supplied individually in plastic wallets.

Surgical Blade No's 6, 9, 10, 10A, 11, 11P, E/11, Sabre E/11  
12, 12D, 13, 14, 15, 15A, 15C, 15T, Sabre D/15, 16, 40 fit  
No's 3, 3L, 5, 7, 9, B3 & B3L handles.

All handles are available in Stainless Steel, individually wrapped, and boxed in units of 10. The easily gripped B3 and B3L are light, adaptable stainless handles, supplied individually in plastic wallets as are the larger 5B Handles.

Surgical Blade No's 18, 19, 20, 21, 22, 22A, 23, 24, 25, 25A, 26, 27, 36 fit No's 4, 4L and 6 handles.

All handles are available in Stainless Steel, individually wrapped, and boxed in units of 10. The larger 6B handle is supplied individually in plastic wallets. *Standard surgical blades in carbon steel (sterile or non-sterile) or in stainless steel (sterile). Boxed in 100's (sterile blades individually foil-wrapped; non-sterile sealed in foil packs of 5). Non-sterile blades also available individually packed in peel-open packets, boxed in 100's. Sabre E/11 and D/15 blades are produced in carbon steel only, sterile, boxed in 100's.*

For use within rigorous Orthopaedic and General surgical procedures. Cygnetic® provides the surgeon with additional strength and stability through the patented combination of its specially designed blade and handle.

Handles are available in Stainless Steel, Non-Sterile, individually sleeved. Blades shapes 10,11,15 and 20 in Stainless Steel, Sterile, individually V.P.I. protected in Metal Foil Packets and boxed in units of 50.

A Post Mortem knife with a choice of two heavy duty carbon steel blades, the pointed PM40 or bull nosed PM40B. Screw type fitting to a stainless steel handle or complete as a convenient disposable scalpel with a protective guard.

Spare blades are individually wrapped and boxed in 10's whilst the scalpels are available in units of 5 and 100.

Further extending the Post Mortem range, the pointed PM60 and bull nosed PM60B blades offer the ease of attaching to a PM8 handle in a traditional manner. Both blades VPI protected in single-peel packets, boxed in 50's. Handles supplied individually in plastic wallets.

Requiring no handles, the standard stitch cutter is in carbon steel, the midi and long versions in stainless steel. The No 3 stainless stitch cutter can be used in conjunction with a No 3 handle for more aseptic techniques.

Available individually foil-wrapped in sterile packs of 100.

Sterilisation is by Gamma Radiation (minimum dosage 25kGy. 2.5 M. Rad) carried out under strict laboratory control in Swan Morton (Services) Ltd's Cobalt 60 continuous irradiation plant.



All medical devices are CE marked in accordance with the Medical Device Directive (93/42/EEC)

"Swann-Morton" and the "ring pattern" logo are registered trademarks of Swann-Morton Limited and related companies. The products shown are manufactured for use by trained health care professionals and can cause injury if used or handled incorrectly. The information and illustration are used in connection with the products shown on this leaflet are for general guidance only and users must seek proper training and instruction in their use before using or handling such products. Swann-Morton Limited exclude all liability for injury to the extent permitted by law.

For cleaning and reprocessing instructions and to view an online video of this procedure, go to <http://www.swann-morton.com/cygnetic>

In stainless steel, sterilised, supplied in single plastic trays in boxes of 10. Fits the Long and Short Swann Major Handles which can be viewed on the website along with the Swann-Major blades.